



# WASH SETS FOR UROLOGY AND ARTHROSCOPY

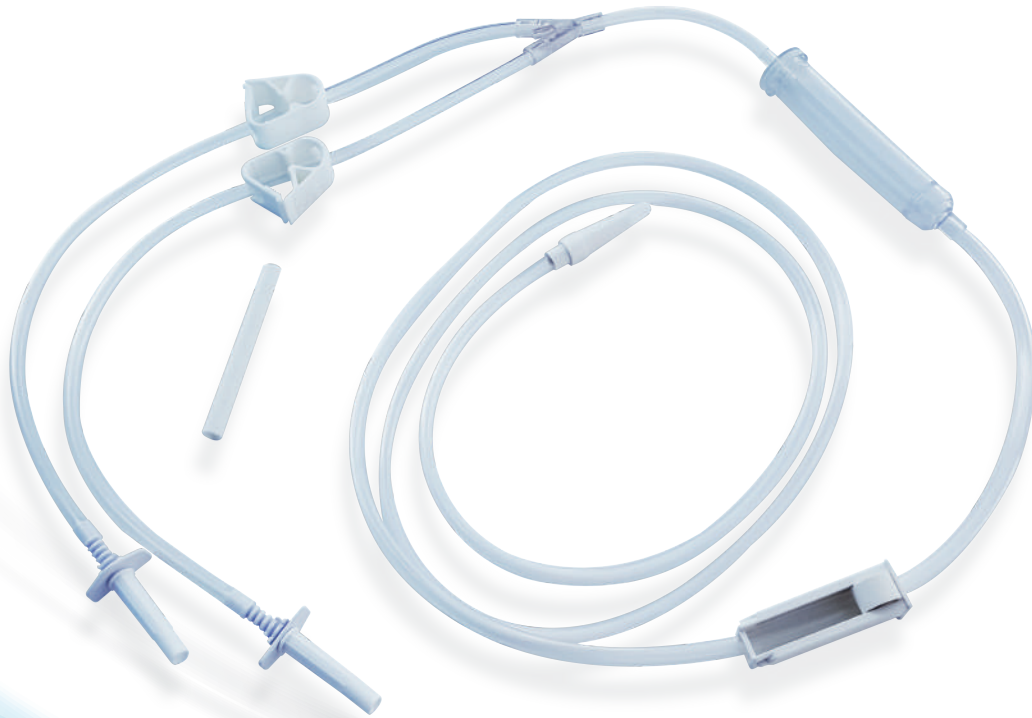

- Tube with wide lumen
- High-flux perforators
- Universal terminal connection




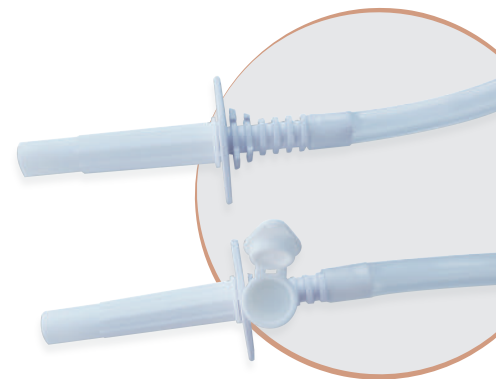


# WASH SETS FOR UROLOGY AND ARTHROSCOPY

Bioservice wash sets, thanks to the wide internal lumen of the tubes, assure high flow rates during urological washing, cystoscopy and arthroscopy procedures. Bioservice manufactures wash sets on specific request of the customer, and has various types of components available, these include:

- 1 ▶ Perforators without filter to be used only for non-rigid containers where emptying does not require external air intake.
- 2 ▶ Perforators with filter which allows for the outflow of fluids from all types of containers: for rigid glass bottles open the filter cap, for non-rigid plastic bags the cap may remain closed.
- 3 ▶ Flow regulators in various dimensions which allow perfect control of the fluid flow rate.
- 4 ▶ Drip chambers in various dimensions which allow for easy checking of the flow rate.
- 5 ▶ Terminal catheter cones with universal connectors fitted with a silicone tube for special connections.
- 6 ▶ Terminal Luer connectors which allow practical and secure attachment to arthroscopes or cystoscopes and easy handling because of the silicone tube into which they are inserted.

The devices are supplied in EtO sterilized packaging.

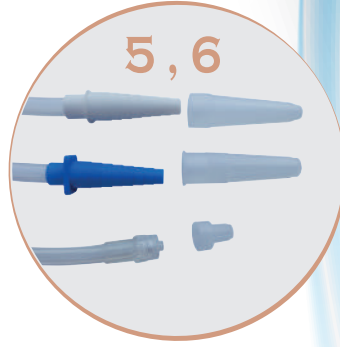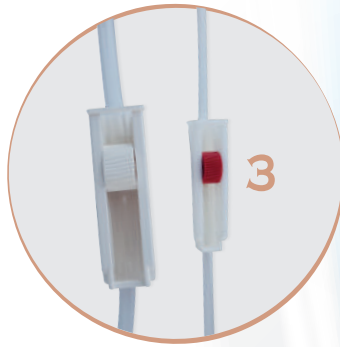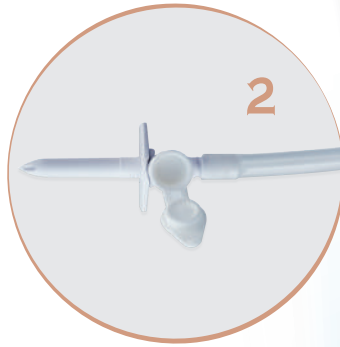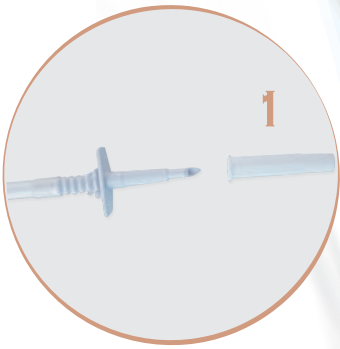

## WASH SETS FOR UROLOGY 150 cmTUBE

CODE	MODEL	COMPOSITION	PACKAGING
H90760E	1-WAY	- 1 perforator without filter with protective cap - drip chamber - roller clamp - terminal catheter cone with protective cap - silicone tube	50 pcs
H90770E	2-WAY	- 2 perforators without filter with protective cap - 2 hose clamps - drip chamber - roller clamp - terminal catheter cone with protective cap - silicone tube	40 pcs
H90780E	3-WAY	- 3 perforators without filter with protective cap - 3 hose clamps - drip chamber - roller clamp - terminal catheter cone with protective cap - silicone tube	30 pcs

## WASH SETS FOR UROLOGY 250 cmTUBE

CODE	MODEL	COMPOSITION	PACKAGING
H90761E	1-WAY	- 1 perforator without filter with protective cap - drip chamber - roller clamp - terminal catheter cone with protective cap	40 pcs
H90771E	2-WAY	- 2 perforators without filter with protective cap - 2 hose clamps - drip chamber - roller clamp - terminal catheter cone with protective cap	40 pcs
H90781E	3-WAY	- 3 perforators without filter with protective cap - 3 hose clamps - drip chamber - roller clamp - terminal catheter cone with protective cap	30 pcs

## WASH SETS FOR ARTHROSCOPY

CODE	MODEL	COMPOSITION	PACKAGING
H90970E	2-WAY	- 2 perforators with filter and protective cap - 2 tube-clamp plates - drip chamber - roller clamp - terminal catheter cone with protective cap - silicone tube	25 pcs
H90980E	2-WAY	- 2 perforators without filter with protective cap - 2 hose clamps - drip chamber - silicone terminal segment with multi-catch clamp - male Luer-lock connector with protective cap	25 pcs