



# CANNULAS AND SETS FOR SURGICAL ASPIRATION

---

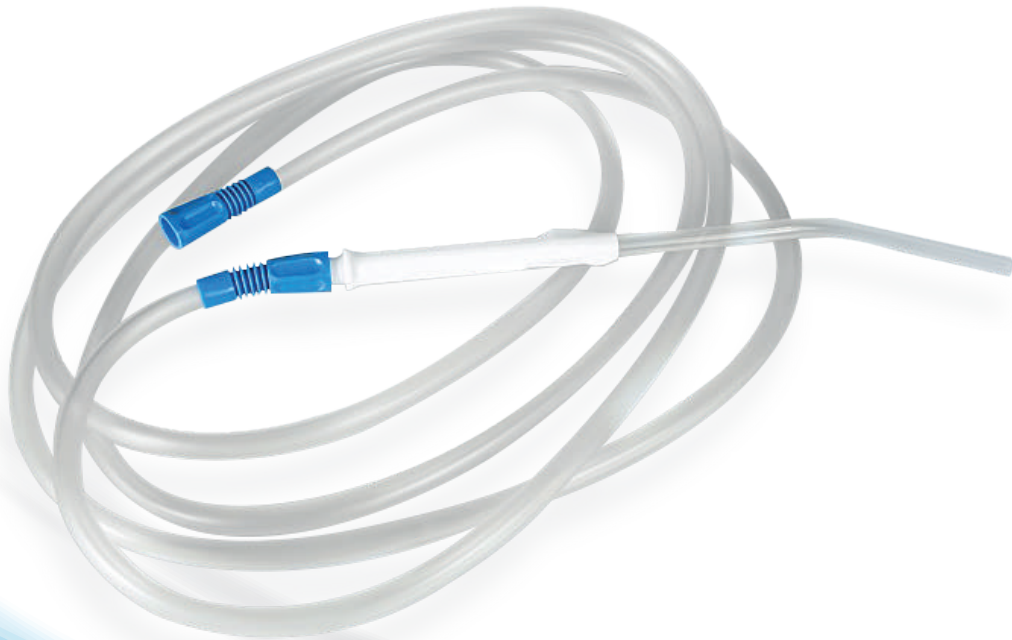

- Cannulas available with and without vacuum control
- Cannulas with atraumatic tip and ergonomic grip
- Soft tube with wide lumen
- Standard connectors

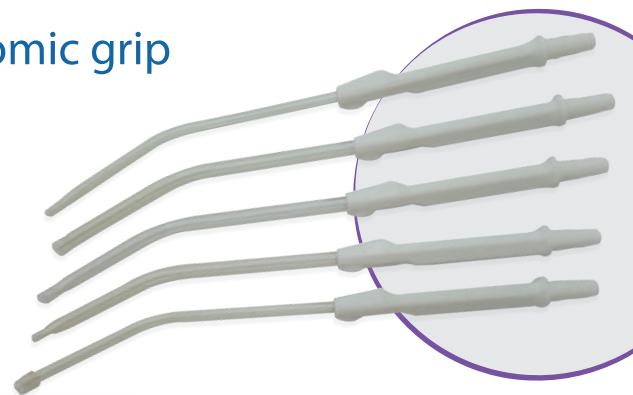

# SURGICAL SUCTION CANNULAS

Yankauer surgical suction cannulas in PVC in various shapes and dimensions.

All the cannula types have been designed in order to satisfy the different market and user requirements.

The cannulas have smooth walls which assure a constant flow of the aspirated fluid and atraumatic tips in various conformations.


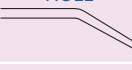





The ergonomic grip makes even the smallest movements easy and safe.

Its lightness makes the device easy to handle.


All versions may also be supplied with vacuum control located on the cannula grip, which allows adjusting aspiration during use.

All the materials used in the production of the suction cannulas have been selected for their technical and qualitative characteristics, and are non-toxic according to the official Pharmacopoeia and are fully biocompatible.

The device is supplied in double packaging sterilized by EtO.

SINGLE CANNULAS						
TIPS AVAILABLE	DIMENSIONS					
	CH12	CH 12 (with vacuum control)	CH20	CH20 (with vacuum control)	CH25	CH25 (with vacuum control)
CURVED DISTAL HOLE 	/	/	H50052E FDC	H50001E VFDC	H50562E FDC	H50511E VFDC
CURVED DISTAL HOLE 	H50351E FDC	H50301E VFDC	/	/	/	/
CURVED DISTAL HOLE 2 lateral holes 	H50358E FDC2FL	H50308E VFDC2FL	/	/	/	/
CURVED DISTAL HOLE 4 lateral holes 	/	/	H50054E FDC4FL	H50003E VFDC4FL	H50564E FDC4FL	H50513E VFDC4FL
FLATTENED CURVED 	/	/	H50055E FDCS	H50004E VFDCS	H50565E FDCS	H50514E VFDCS
FLATTENED CURVED 2 lateral holes 	/	/	H50056E FDCS2FL	H50005E VFDCS2FL	H50566E FDCS2FL	H50515E VFDCS2FL
SOFT CURVED GRID 	/	/	H50057E CCM	H50006E VCCM	H50567E CCM	H50516E VCCM

All the cannulas are available in multi-packs of 50 pieces.

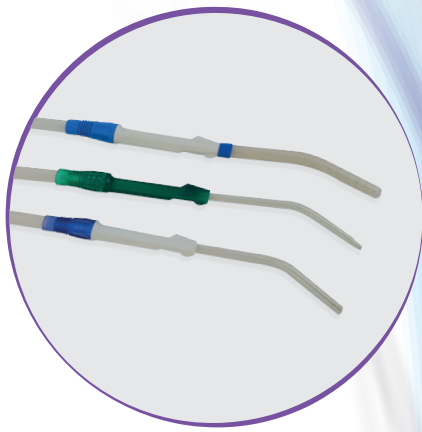

# SURGICAL ASPIRATION SET

The set is composed of a cannula connected to the suction tube.

The tubes are extremely soft and have a wide lumen to guarantee excellent flow of the aspirated fluid.



The tube has been designed in such a way as to prevent lumen constrictions during the aspiration phases.

The tube is available in different standard lengths or customisable on the basis of the customer's requirements.

The tube is equipped with female connectors whose bellows geometry allows easy handling and bending over 90° without reducing the lumen.

The connectors guarantee perfect seal.

The device is supplied in double packaging sterilized by EtO.

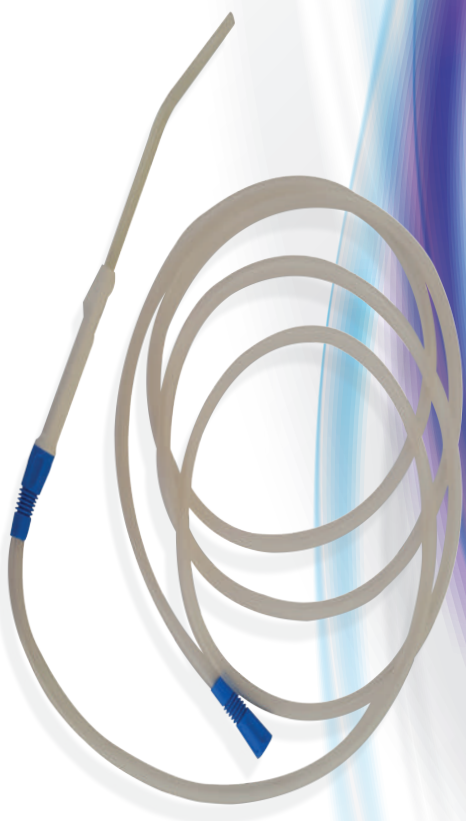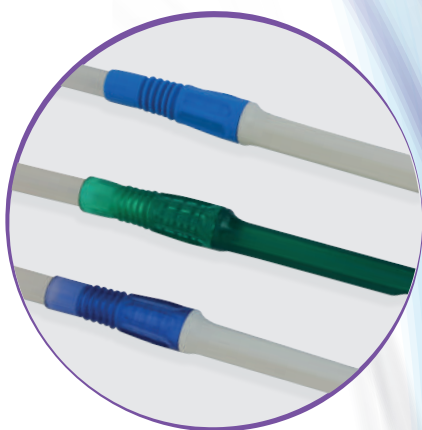

CH12 CANNULA+ CH25 TUBE				
200 cm TUBE			300 cm TUBE	
TIPS AVAILABLE	CH12 CANNULA	CH12 CANNULA (with vacuum control)	CH12 CANNULA	CH12 CANNULA (with vacuum control)
CURVED DISTAL HOLE	H54211 E	H54020 E	H54215 E	H54024 E
CURVED DISTAL HOLE 2 lateral holes	H54323 E	H54132 E	H54327 E	H54136 E

CH20 CANNULA+ CH25 TUBE				
200 cm TUBE			300 cm TUBE	
TIPS AVAILABLE	CH20 CANNULA	CH20 CANNULA (with vacuum control)	CH20 CANNULA	CH20 CANNULA (with vacuum control)
CURVED DISTAL HOLE	H51211 E	H51020 E	H51215 E	H51024 E
CURVED DISTAL HOLE 4 lateral holes	H51243 E	H51052 E	H51247 E	H51056 E
FLATTENED CURVED	H51259 E	H51068 E	H51263 E	H51072 E
FLATTENED CURVED 2 lateral holes	H51275 E	H51084 E	H51279 E	H51088 E
SOFT CURVED GRID	H51291 E	H51100 E	H51295 E	H51104 E

CH25 CANNULA+ CH30 TUBE				
200 cm TUBE			300 cm TUBE	
TIPS AVAILABLE	CH25 CANNULA	CH25 CANNULA (with vacuum control)	CH25 CANNULA	CH25 CANNULA (with vacuum control)
CURVED DISTAL HOLE	H52595 E	H52403 E	H52599 E	H52407 E
CURVED DISTAL HOLE 4 lateral holes	H52627 E	H52435 E	H52631 E	H52439 E
FLATTENED CURVED	H52643 E	H52451 E	H52647 E	H52455 E
FLATTENED CURVED 2 lateral holes	H52659 E	H52467 E	H52663 E	H52471 E
SOFT CURVED GRID	H52675 E	H52483 E	H52679 E	H52487 E

All the surgical aspiration sets are available in multi-packs of 25 pieces.

# CANNULAS AND SETS FOR SURGICAL ASPIRATION BIOPOOL



## SUCTION CANNULAS

"Biopool" was conceptualized and created out of the need to perform aspiration at locations difficult to reach when using standard cannulas.

"Biopool" is made up of an internal cannula and flexible perforated external cannula with an atraumatic tip. By sliding the external cannula over the internal one, the desired extension and aspiration flow adjustment is obtained. The holes in the external cannula guarantee high aspirating capacity. The ergonomic grip makes even the smallest movements easy and safe. Its lightness makes the devices easy to handle. All the versions may be supplied with vacuum control.

All the materials used have been selected for their technical and qualitative characteristics, and are non-toxic according to the Official Pharmacopoeia and the fully biocompatible.

The device is supplied in double packaging sterilized by EtO.

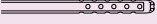SINGLE CANNULAS		
TIPS AVAILABLE	DIMENSIONS	
	CH20	CH20 (with vacuum control)
STRAIGHT 	H50062E POOL A	H50011E V POOL A
CURVED 10° 	H50064E C POOL A	H50013E VC POOL A

All the cannulas are available in multi-packs of 50 pieces.

## ASPIRATION SET

The set is composed of a PVC cannulas connected to the suction tube. The tubes are extremely soft and have a wide lumen to guarantee excellent flow of the aspirated fluid. The tube has been designed in such a way as to prevent lumen constrictions during the aspiration phases. The tube is available in different standard lengths or customisable on the basis of the customer's requirements.

The device is supplied in double packaging sterilized by EtO.

CANNULAS + CH25 TUBE (LENGTH 200CM)		
TIPS AVAILABLE	CANNULA DIMENSIONS	
	CH20 CANNULA	CH20 CANNULA (with vacuum control)
STRAIGHT	H51371E	H51179E
CURVED 10°	H51381E	H51189E

CANNULAS + CH25 TUBE (LENGTH 300CM)		
TIPS AVAILABLE	CANNULA DIMENSIONS	
	CH20 CANNULA	CH20 CANNULA (with vacuum control)
STRAIGHT	H51375E	H51183E
CURVED 10°	H51382E	H51190E

CANNULAS + CH30 TUBE (LENGTH 200CM)		
TIPS AVAILABLE	CANNULA DIMENSIONS	
	CH20 CANNULA	CH20 CANNULA (with vacuum control)
STRAIGHT	H51755E	H51563E
CURVED 10°	H51765E	H51573E

CANNULAS + CH30 TUBE (LENGTH 300CM)		
TIPS AVAILABLE	CANNULA DIMENSIONS	
	CH20 CANNULA	CH20 CANNULA (with vacuum control)
STRAIGHT	H51759E	H51567E
CURVED 10°	H51766E	H51574E

All the surgical aspiration sets are available in multi-packs of 25 pieces.